

PRECISE BIOPHARMA PVT. LTD.

Product: Semaglutide Injection 2 mg/3 ml					
Product Name: CELTUIDE			Customer Name: TaxaneHealth Care (Cambodia) Co., Ltd.		
Open Size: L-297 x H-420 mm		Size After Fold: NA		Style: Pack Insert (IFU)	
Market: Export		Country: Cambodia		Colour: Black & White (CMYK)	
Substrate: NA					
Created on: 21-11-2025		Last Revised on: 24-11-2025		Item Code: KH-TID2-SI-IFU-00	
Country Code: NA					
Manufacturing Location: Sovereign Pharma Private Limited					
Checked by	Artist	RA	QA	Plant QA	Country QA

297.00 mm

INSTRUCTION FOR USE

**Semaglutide Injection 2mg/3ml**

for subcutaneous use

**CELTUIDE**

**0.25 mg or 0.5 mg doses**

(pen delivers doses in 0.25 mg or 0.5 mg increments only)


- Read these instructions carefully before using your Semaglutide Injection 2mg/3ml pen.
- Do not use your pen without proper training from your healthcare provider. Make sure that you know how to give yourself an injection with the pen before you start your treatment.
- Do not share your Semaglutide Injection 2mg/3ml pen with other people, even if the needle has been changed. You may give other people a serious infection, or get a serious infection from them.

**⚠ If you are blind or have poor eyesight and cannot read the dose counter on the pen, do not use this pen without help.** Get help from a person with good eyesight who is trained to use the Semaglutide Injection 2mg/3ml pen.

- Start by checking your pen to make sure that it contains Semaglutide Injection, and then look at the pictures below to get to know the different parts of your pen and needle.
- Your pen is a prefilled, single-patient-use; dial-a-dose pen. It contains 2 mg of semaglutide, and you can select doses of 0.25 mg or 0.5 mg. Each prefilled pen contains 4 doses of 0.25 mg and 2 doses of 0.5 mg or contains 4 doses of 0.5 mg.
- Your pen is made to be used with disposable needles up to a length of 8 mm.
- 32G 4mm disposable needles are included with your Semaglutide Injection 2mg/3ml pen.
- Always use a new needle for each injection.
- Supplies you will need to give your Semaglutide injection:
  - Pen
  - A new needle


See "Disposing of used Semaglutide Injection pens and needles" at the end of these instructions.

**Semaglutide Pen and Needles (Example)**

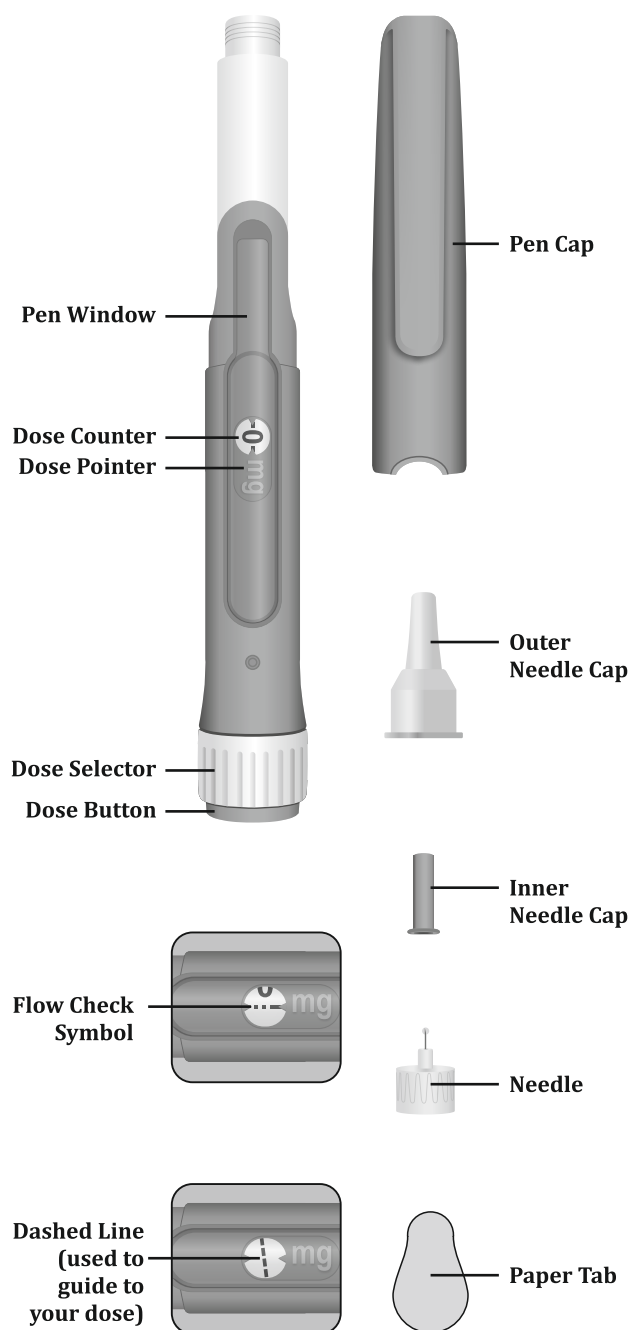


**Step 1. Prepare your pen with a new needle**


- Wash your hands with soap and water.
- Check the name and label of your pen, to make sure that it contains Semaglutide Injection 2mg/3ml. This is especially important if you take more than 1 type of medicine.
- Pull off the pen cap.

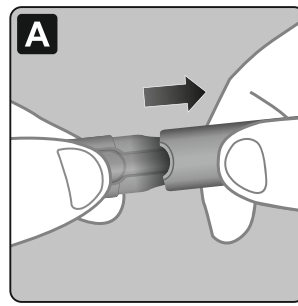


- Check that the Semaglutide Injection 2mg/3ml medicine in your pen is clear and colorless. Look through the pen window. If Semaglutide Injection 2mg/3ml looks cloudy or contains particles, do not use the pen.

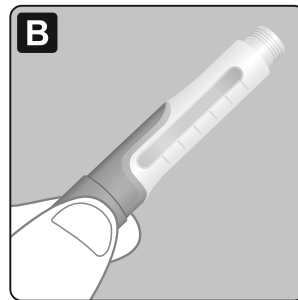


- Take a new needle, and tear off the paper tab. Do not attach a new needle to your pen until you are ready to give your injection.

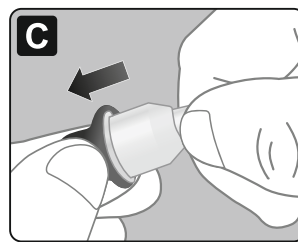


- Push the needle straight onto the pen. Turn until it is on tight.

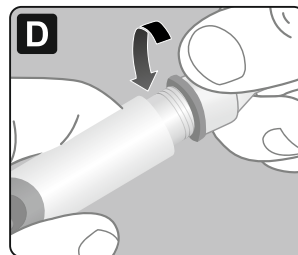

- Pull off the outer needle cap. Do not throw it away.

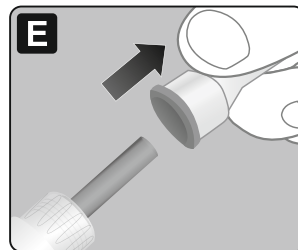

- Pull off the inner needle cap and throw it away. A drop of Semaglutide Injection may appear at the needle tip. This is normal, but you must still check the Semaglutide Injection flow if you use a new pen for the first time.

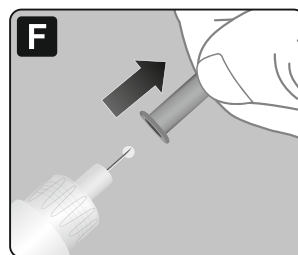


Always use a new needle for each injection. This will reduce the risk of contamination, infection, leakage of Semaglutide Injection 2mg/3ml, and blocked needles leading to the wrong dose.


**⚠ Do not reuse or share your needles with other people. You may give other people a serious infection, or get a serious infection from them. Never use a bent or damaged needle.** Do not attach a new needle to your pen until you are ready to take your injection.

**Step 2. Check the Semaglutide Injection flow with each new pen**

- Check the Semaglutide Injection 2mg/3ml flow before the first injection with each new pen only.
- If your Semaglutide Injection 2mg/3ml pen is already in use, go to Step 3 "Select your dose".
- Turn the dose selector until the dose counter shows the flow check symbol (---)

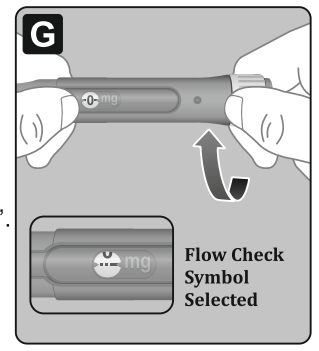


- Hold the pen with the needle pointing up. Press and hold in the dose button until the dose counter shows 0. The 0 must line up with the dose pointer. A drop of Semaglutide Injection 2mg/3ml will appear at the needle tip.
- If no drop appears, repeat Step 2 above as shown in Figure G and Figure H up to 6 times. If there is still no drop, change the needle and repeat Step 2 as shown in Figure G and Figure H 1 more time.
- Do not use the pen if a drop of Semaglutide Injection 2mg/3ml still does not appear.

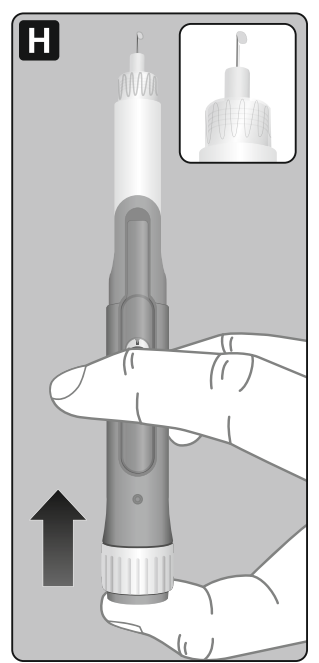

**⚠ Always make sure that a drop appears at the needle tip before you use a new pen for the first time.** This makes sure that Semaglutide Injection 2mg/3ml flows. If no drop appears, you will not inject any Semaglutide Injection 2mg/3ml, even though the dose counter may move. This may mean that there is a blocked or damaged needle. A small drop may remain at the needle tip, but it will not be injected. Only check the Semaglutide Injection 2mg/3ml flow before your first injection with each new pen.

**Step 3. Select your dose**


- Turn the dose selector until the dose counter stops and shows your dose (0.25 mg or 0.5 mg). The dashed line in the dose counter (|)

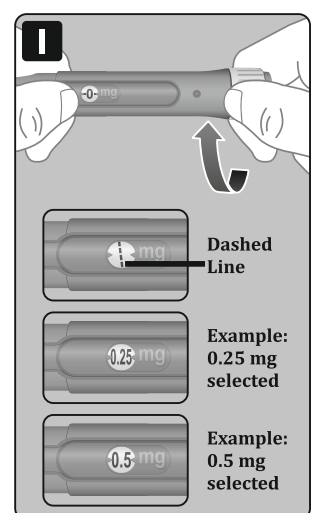

**⚠ Always use the dose counter and the dose pointer to see how many mg you select.** You will hear a "click" every time you turn the dose selector. Do not set the dose by counting the number of clicks you hear. Only doses of 0.25 mg or 0.5 mg can be selected with the dose selector. The selected dose must line up exactly with the dose pointer to make sure that you get a correct dose. The dose selector changes the dose. Only the dose counter and dose pointer will show how many mg you select for each dose. You can only select 0.25 mg or 0.5 mg for each dose. When your pen contains less than 0.25 mg or 0.5 mg, the dose counter stops before 0.25 mg or 0.5 mg is shown. The dose selector clicks differently when turned forward or backward. Do not count the pen clicks.






### How much Semaglutide Injection 2mg/3ml is left?

- To see how much Semaglutide Injection 2mg/3ml is left in your pen, use the dose counter: Turn the dose selector until the **dose counter stops**.
- If it shows 0.5, **at least 0.5 mg** is left in your pen. If the **dose counter stops before 0.5 mg**, there is not enough Semaglutide Injection 2mg/3ml left for a full dose of 0.5 mg.
- If it shows 0.25, **at least 0.25 mg** is left in your pen. If the **dose counter stops before 0.25 mg**, there is not enough Semaglutide Injection 2mg/3ml left for a full dose of 0.25 mg.
- If there is not enough Semaglutide Injection 2mg/3ml left in your pen for a full dose, do not use it. Use a new Semaglutide Injection 2mg/3ml pen.

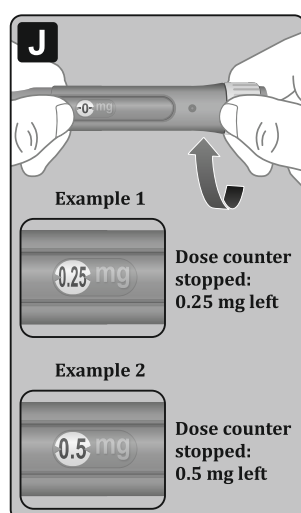

### Step 4. Inject your dose

- Choose your injection site and wipe the skin with an alcohol swab. Let the injection site dry before you inject your dose. (See Figure K)
- Insert the needle into your skin as your healthcare provider has shown you. **Make sure you can see the dose counter.** Do not cover it with your fingers. This could stop the injection.
- Press and hold down the dose button until the dose counter shows 0. The 0 must line up with the dose pointer. You may then hear or feel a click.
- Keep the needle in your skin after the dose counter has returned to 0 and count slowly to 6. If the needle is removed earlier, you may see a stream of Semaglutide Injection 2mg/3ml coming from the needle tip. If this happens, the full dose will not be delivered.
- Remove the needle from your skin. If blood appears at the injection site, press lightly with a gauze pad or cotton ball. Do not rub the area.

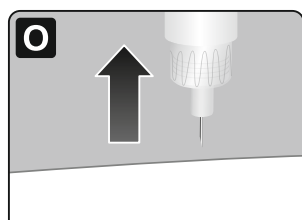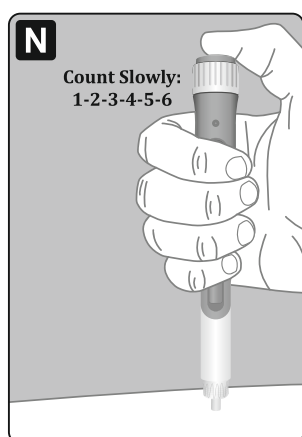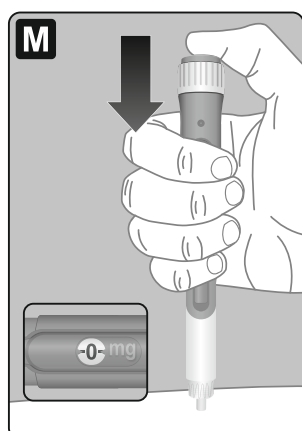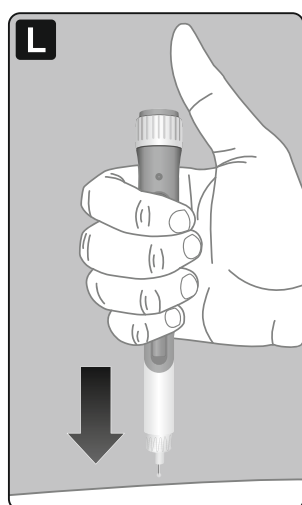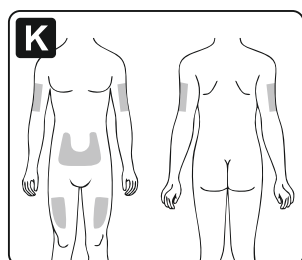

- Always watch the dose counter to make sure you have injected your complete dose.** Hold the dose button down until the dose counter shows 0.

#### How to identify a blocked or damaged needle?

- If 0 does not appear in the dose counter after continuously pressing the dose button, you may have used a blocked or damaged needle.
- If this happens you have **not received any** Semaglutide Injection 2mg/3ml even though the dose counter has moved from the original dose that you have set.


#### How to handle a blocked needle?

- Change the needle as described in Step 5, and repeat all steps starting with Step 1: "Prepare your pen with a new needle".
- Never touch the dose counter when you inject.** This can stop the injection.


You may see a drop of Semaglutide Injection 2mg/3ml at the needle tip after injecting. This is normal and does not affect your dose.

### Step 5. After your injection


- Carefully remove the needle from the pen. Do not put the needle caps back on the needle to avoid needle sticks.

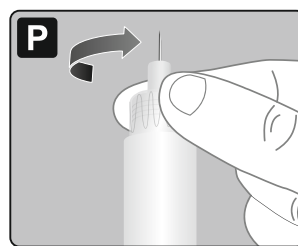


- Place the needle in a sharps disposal container right away to reduce the risk of needle sticks. See "Disposing of used Semaglutide Injection 2mg/3ml pens and needles" below for more information about how to dispose of used pens and needles the right way.

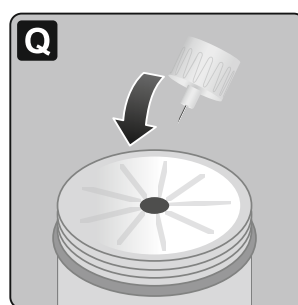

- Put the pen cap on your pen after each use to protect Semaglutide Injection 2mg/3ml from light.

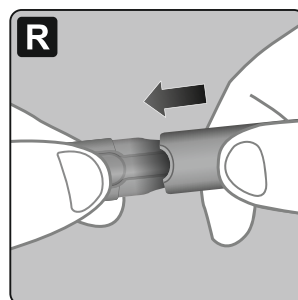

- If you do not have a sharps disposal container, follow a 1-handed needle recapping method. Carefully slip the needle into the outer needle cap. Dispose of the needle in a sharps disposal container as soon as possible.

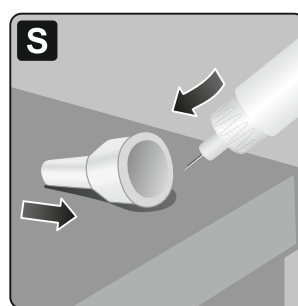

- Never try to put the inner needle cap back on the needle.** You may stick yourself with the needle. **Always remove the needle from your pen.** This will reduce the risk of contamination, infection, leakage of Semaglutide Injection 2mg/3ml, and blocked needles leading to the wrong dose. If the needle is blocked, you will not inject any Semaglutide Injection 2mg/3ml. **Always dispose of the needle after each injection.**

#### Disposing of used Semaglutide Injection 2mg/3ml pens and needles:

- Put your used Semaglutide Injection 2mg/3ml pen and needle in a FDA-cleared sharps disposal container right away after use.
- If you do not have a FDA-cleared sharps disposal container, you may use a household container that is:
  - made of a heavy-duty plastic
  - can be closed with a tight-fitting, puncture-resistant lid, without sharps being able to come out
  - upright and stable during use
  - leak-resistant
  - properly labeled to warn of hazardous waste inside the container
- When your sharps disposal container is almost full, you will need to follow your community guidelines for the right way to dispose of your sharps disposal container. There may be state or local laws about how you should throw away used needles and syringes.

- Do not dispose of your used sharps disposal container in your household trash unless your community guidelines permit this. Do not recycle your used sharps disposal container.
- Safely dispose of Semaglutide Injection 2mg/3ml that is out of date or no longer needed.

### Important

- Caregivers must be **very careful when handling used needles** to prevent accidental needle stick injuries and prevent passing (transmission) of infection.
- Never use a syringe to withdraw Semaglutide Injection 2mg/3ml from your pen.
- Always carry an extra pen and new needles** with you, in case of loss or damage.
- Always keep your pen and needles out of reach of others**, especially children.
- Always keep your pen with you.** Do not leave it in a car or other place where it can get too hot or too cold.

### Caring for your pen

- Do not drop your pen** or knock it against hard surfaces. If you drop it or suspect a problem, attach a new needle and check the Semaglutide Injection 2mg/3ml flow before you inject.
- Do not try to repair your pen** or pull it apart.
- Do not expose your pen to dust, dirt or liquid.**
- Do not wash, soak, or lubricate your pen.** If necessary, clean it with mild detergent on a moistened cloth.

### How should I store my Semaglutide Injection 2mg/3ml pen?

- Store your new, unused Semaglutide Injection 2mg/3ml pens in the refrigerator between 2°C to 8°C.
- Store your pen in use** for 56 days below 30°C or in a refrigerator between 2°C to 8°C.
- The Semaglutide Injection 2mg/3ml pen you are using should be disposed of (thrown away) after 56 days, even if it still has Semaglutide Injection 2mg/3ml left in it. Write the disposal date on your calendar.
- Do not freeze Semaglutide Injection 2mg/3ml. Do not use Semaglutide Injection 2mg/3ml if it has been frozen.
- Unused Semaglutide Injection 2mg/3ml pens may be used until the expiration date ("EXP") printed on the label, if kept in the refrigerator.
- When stored in the refrigerator, do not store Semaglutide Injection 2mg/3ml pens directly next to the cooling element.
- Keep Semaglutide Injection 2mg/3ml away from heat and out of the light.
- Keep the pen cap on when not in use.
- Keep Semaglutide Injection 2mg/3ml and all medicines out of the reach of children.

#### Marketed By :

**Taxane**  
TAXANE HEALTH CARE  
(CAMBODIA) CO., LTD.

**TaxaneHealth Care (Cambodia) Co., Ltd**

PT-1258 st prek tasek borey aphivath town (ocic), Daeum Kor, Chrouy Changvar, Phnom Penh, Cambodia.  
+855 96 309 6850

#### Manufactured At :

**Sovereign Pharma Pvt. Ltd.**

Survey No. 46/1-4, Kadaiya Village, Daman-396 210, India.  
Mfg. Lic. No. : DD/L/868

#### Developed By :

Precise Biopharma Pvt. Ltd., India.

04-T02-SI-RU-00